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
Electric toothbrush

Instructions

Tips for using

- Please note that due to its effectiveness, the OrisCura is not recommended for use by children under12 years of age!
- The Oris Cura electric toothbrush is the modern toothbrush of the new era, powered by an elongated main brush.
- The cleaning effect of the Oris Cura electric toothbrush is based on a triple action. The soft bristles of the brush head and the rotating movement of the bristles clean your teeth in a gentle way. In addition, the elongatedsection of the brush head guides a mixture of saliva, water and toothpaste into the interdental spaces and cleans them effectively and prevents damage to the gums.
- To operate the Oris Cura electric toothbrush according to the invention, the handle may include one or more push buttons that control the rotation of the brush section and the revolutions per minute. The buttons on the handle are adapted to control either the upper or lower jaw. The direction of rotation is then automatically the correct one, i.e. brushing is done from the gums to the teeth. The handle may, for example, include two push buttons, where one push button adjusts the direction of rotation for the upper jaw and the other button adjusts the direction of rotation for the lower jaw so that brushing always takes place from the gums to the teeth.

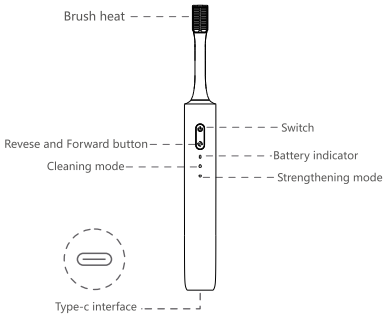Product Material	Main ABS material+DuPont brush wire
Vibration revolution	300 rpm
Waterproof grade	IPX7 level
Charging time	3-4 hours
Usage duration	Normally usable for 3 months

* The usage time is calculated based on 3 minutes per session and 2 times per day

Function Introduction

- On/Off
Press the switch button to start the machine, then press the power button to switch gears, and press any gear long to shut down. The default working time is 3 minutes, and the machine will shut down automatically after 3 minutes. The machine has a gear memory function, which defaults to the last used mode when it is turned on next time.
- Two major modes
Clean-Clean mode: for daily cleaning use.
Strong-strengthening mode: suitable for stubborn tooth stains.
- Forward and reverse mode
Press the forward and reverse button to switch the rotation direction.
- Battery display
When the power is low, the red light flashes, the red light breathes when charging, and the green light is steady on after full. When the power does not meet the startup voltage, the red light is steady on for 3s, and the shutdown does not work.

- We recommend that you start getting used to the Oris cura with the "Sensitive" program (100 revolutions per minute). The microbubbles and the elongated brush section may be itchy at first, but you get used to it after a few brushings.
- The electric toothbrush is moved in the same way as a regular toothbrush, but slower. The brush is placed on the gum line, even if you touch the gum line, there is no danger at all because the direction of rotation is set correctly. The brush does NOT need to be moved in a circular oval motion as you have to do with other toothbrushes.
- According to the design, the electric toothbrush has a guard on the back of the toothbrush head to protect the inner cheek. The electric toothbrush has a function where the brush itself either rotates downwards or upwards depending on where the toothbrush is placed, on the lower jaw or the upper jaw.

- The elongated movement of the brushes provides an effective and gentle cleaning result, even in hard-to-reach places and along the gum line. The brush head has a narrow neck to reach further and remove plaque.

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